


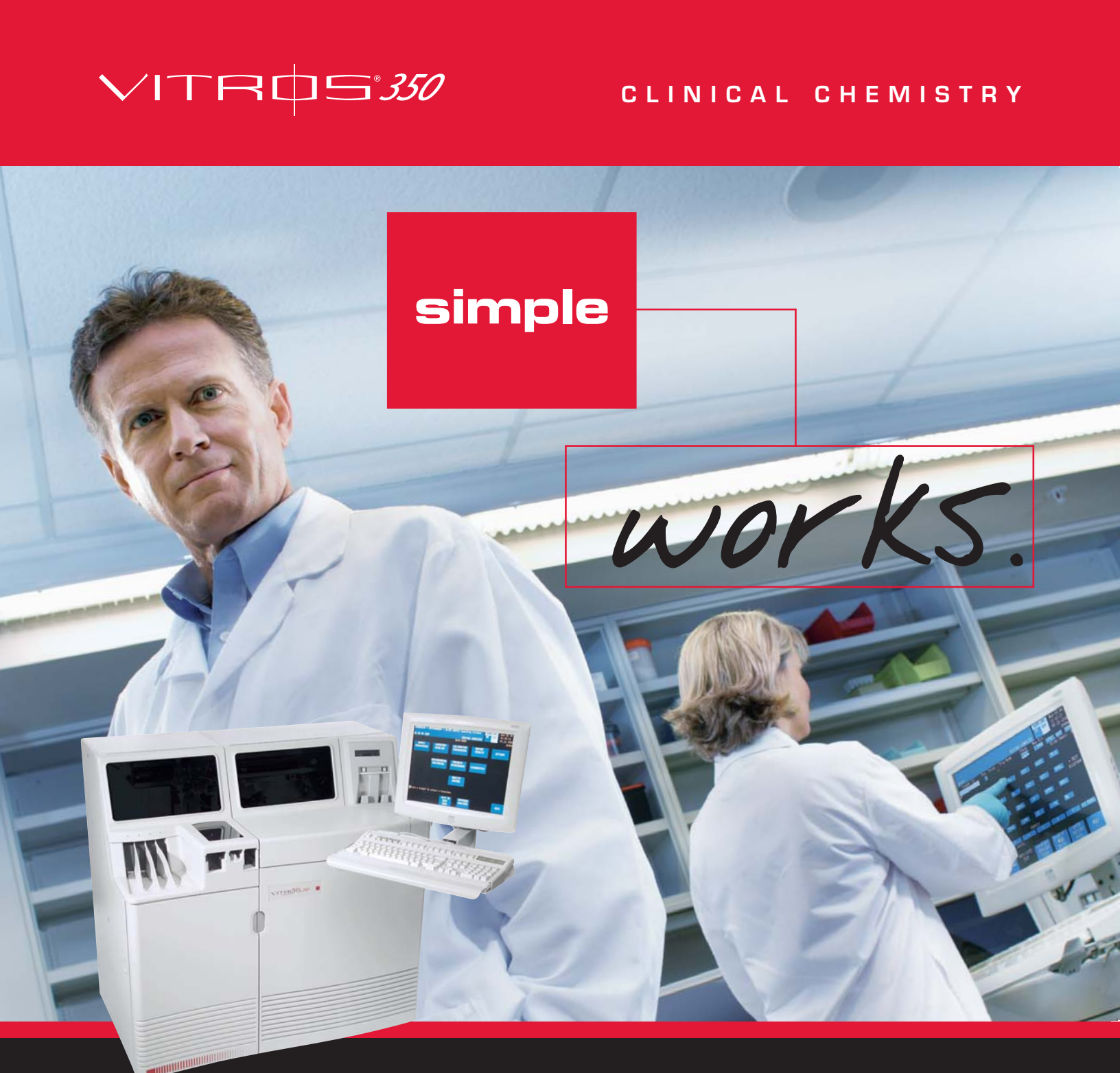

simple

works.

The new VITROS® 350 Chemistry System is designed to be efficient, dependable, and — above all — easy to use. **FACT:** In a recent survey of 535 laboratory professionals, “simplified operation that reduces error” was the most frequently chosen (73%) and valued criteria in selecting a new chemistry system.* And customers rate VITROS Chemistry Systems as the easiest to use and most reliable.** Whether you use it as a primary, STAT, or backup system, the VITROS 350 Chemistry System will get the job done right the first time. **IT'S THAT SIMPLE.** For more information, contact your Ortho-Clinical Diagnostics Representative or visit our Web site at www.orthoclinical.com.

* CaseBauer™ Chemistry and Immunoassay Analyzers, North America and Europe 2002

** EAC U.S. Hospital Clinical Chemistry and Immunoassay Survey 2003 (VITROS 950 and 250 Systems) and Frost & Sullivan European Clinical Laboratories Analysers Market 2004



Ortho-Clinical Diagnostics

a Johnson & Johnson company

VITROS® 350 Chemistry System Specifications

Labor Saving Ease of Use

- The VITROS 350 doesn't require highly-trained operators—simple procedures allow around-the-clock availability to all members of the lab staff
- Easy to train and cross train operators
- Walk-away operation
- Quality control procedures are required just once each day and calibration intervals up to 6 months
- Process monitoring

Operator Interface

- Ergonomically designed, flat low-glare LCD touch screen monitor
- Keyboard platform and support arm
- Flexible, customized positioning
- On-board documentation and help function

Versatile and Efficient

- Up to 300 results per hour
- On-board capacity for all 45 assays
- Primary tube sampling
- On-board dilution
- Optional positive patient identification via barcode scanner
- Flexible processing for manual, stat or automated testing

Secure Results


- Unique sample processing capabilities like disposable tips, bubble and clot detection assure confidence in results
- Self-monitoring alerts operators to prompt action, avoiding testing delays

Environment Friendly

- No plumbing, drains, vents, deionized water is required; all “waste” is contained in used test slides which are disposed daily

Cost-effective VITROS MicroSlide™ Technology

Delivers low cost per reportable result and high reagent efficiency without the costs, maintenance, preparation, carryover, and interference associated with traditional water-based and indirect ISE systems.

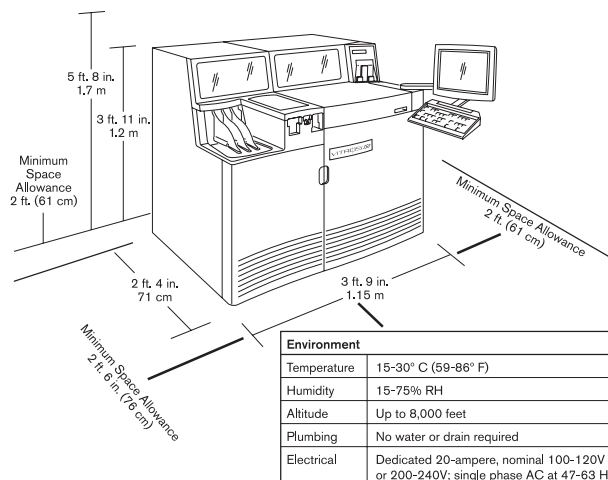

- **Measurement Principles:** Potentiometric (direct ISE), Colormetric/Rate, Immuno-rate
- **On-Board Test Capacity:** 3600 tests
- **Sample Types:** Serum, Plasma, Urine, CSF
- **Sample Volume:** 5-11 µL
- **Sample Capacity:** 40 samples (4 trays of 10 samples)
- **Time for Single Results:** ~ 2.5 to 8 minutes
 - Potentiometric: ~ 2.5 minutes
 - Colormetric: ~ 5 minutes
 - Immuno-rate: ~ 8 minutes
- **Sample Containers:** 10mL, 7mL, 5mL, 2-4mL Collection Tubes, 1mL Microtube, Microsample cup and .5mL and 2mL
- **Sample Management:** Clot Detection, Bubble Detection, Liquid Level-Sensing (Pressure Transducer), Short-Sample Detection
- **Automated Technology (optional):** Includes automatic tip loading and deep tube sampling
- **Calibration:** Extended; multiple lots supported

VITROS MicroSlide Tests

Acetaminophen	Cholesterol/HDL	Lactate	Total Bilirubin
Alanine Aminotransferase	Chloride	Lactate Dehydrogenase	Total Protein
Albumin	Cholinesterase	Lipase	Triglyceride
Alcohol	Creatinine	Lithium	Unconjugated/Conjugated Bilirubin
Alkaline Phosphatase	Creatine Kinase	Magnesium	Urea Nitrogen
Ammonia	Creatine Kinase MB	Phenobarbital	Uric Acid
Amylase	CSF Protein	Phenytoin	Urine Protein
Aspartate Aminotransferase	Digoxin	Phosphorus	
C-Reactive Protein	Gamma Glutamyl Transferase	Potassium	
Calcium	Glucose	Salicylate	
Carbamazepine	dHDL*	Sodium	
Carbon Dioxide	Iron/Total Iron Binding Capacity	Theophylline	

*under development

For more information contact your Ortho-Clinical Diagnostics representative or visit our website at www.orthoclinical.com.